

# iOSBrainLinkProSDKSpecialEdition

---

## V1.3\_Release

Date: 1116, 2020

Author: Liang Fang

SDK Version: 1.3



### **Update record:**

1. V1.3 fix the bug of failed callback in iOS13.5 system
2. V1.2 supports multiple or single connections
3. V1.0 supports single connection

# Content

---

iOSBrainLinkProSDKSpecialEditon V1.3_Release .....	1
iOSBrainLinkProSDKSpecialEditon Development Guide.....	5
Introduction.....	5
Your First Project: iOSBrainLinkProSDKSpecialEditon V1.3_Release	
Demo .....	6
iOSBrainLinkProSDKSpecialEditon V1.3 API Reference .....	12
HZLBlueData Reference .....	12
Blue4ManagerReference .....	14

# iOSBrainLinkProSDKSpecialEdition Development

---

## Guide

### Introduction

This guide will teach you how to use HZLBlueTooth SDK for iOS to write iOS applications that can acquire brainwave data from MacroTelligence 's Hardware ( BrainLink Pro & BrainLink Lite ) . This will enable your iOS apps to receive and use brainwave data such as BLEMIND and BLEGRAVITY acquired via Bluetooth, MacroTelligence 's Hardware and File source encapsulated as HZLBlueTooth. HZLBlueTooth SDK for iOS supports upgrading Hardware

Function:

Receive brainwave data. One or more Bluetooth devices can be connected at the same time.

**Files included:**

- API Reference (this document)
- SDK static library and headers
- iOSBrainLinkProSDKSpecialEdition V1.3\_Release.a
- HZLBlueData.h
- Blue4Manager.h
- iOSBrainLinkProSDKSpecialEdition V1.3\_Release Demo

**Supported devices:**

- Data format with power
  - BrainLink\_Pro
  - Jii
- Data format without power
  - BrainLink\_Lite
  - Mind Link

**iOS VerrSION :**


- iOS 9.0 +

# Your First Project: iOSBrainLinkProSDKSpecialEdition


## V1.3\_Release Demo

### Step 1:

1.1 Import the IOS framework libraries CoreBluetooth.framework in the Build Phases of TARGETS in the Xcode project:

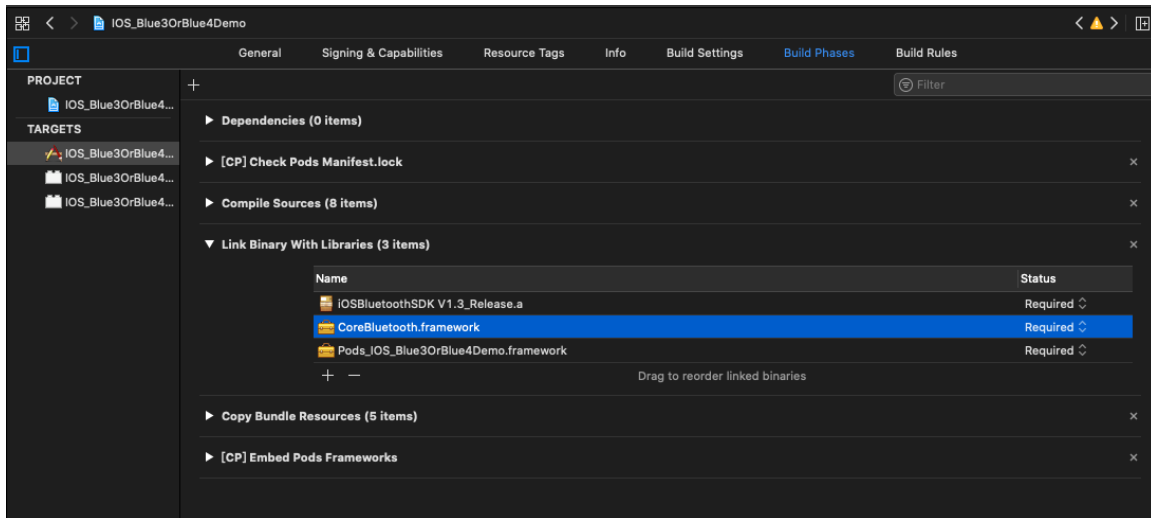


Add the Bluetooth permission privacy permission privacy In the Info.plist: ( ios13 needs to add the Bluetooth permission privacy - Bluetooth always usage description and privacy - Bluetooth peripheral usage description )

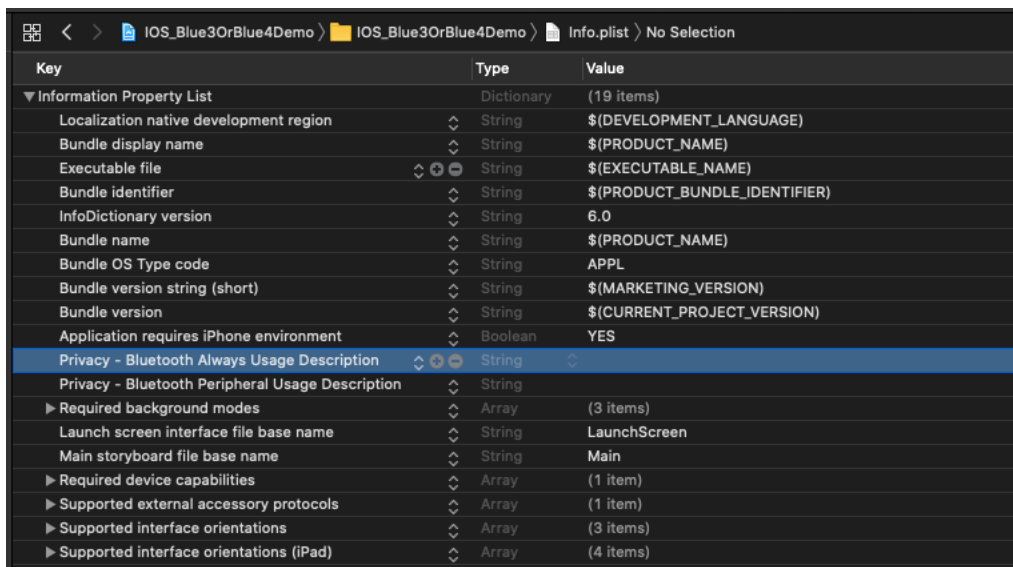



1.2 If you want Bluetooth to work in the background, please set it as follows ,  
Don't set it if you don't need it:

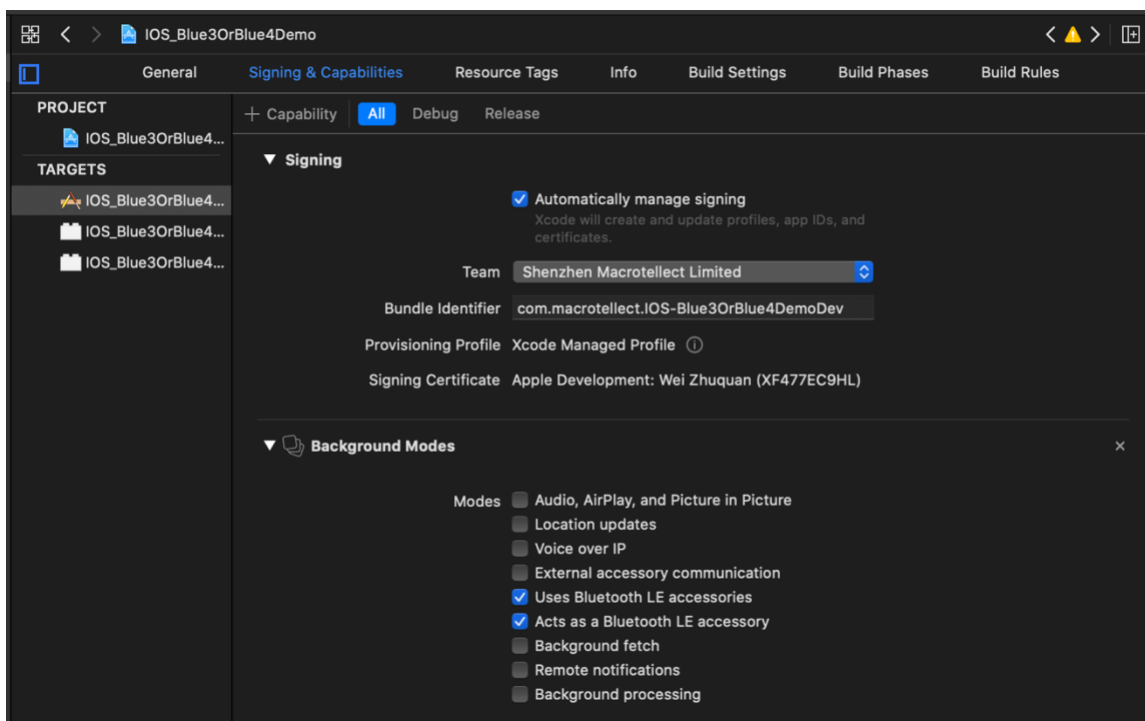

## Step 2:

Import header file

```
#import "HZLBlueData.h"
```

```
#import "Blue4Manager.h"
```

Function one : Receive data

```
// Bluetooth multi-connection
```

```
// blueNames: connectable device
```

```
NSArray *blueNames = @[@"BrainLink",@"BrainLink_Pro",@"jii@jii-***"];
```

```
[Blue4Manager logEnable:YES];
```

```
[[Blue4Manager sharedInstance] configureBlueNames:blueNames ableDeviceSum:7];
```

```
// Connect bluetooth callback successfully
```

```
__weak FactoryViewController *weakSelf = self;
```

```
[Blue4Manager sharedInstance].blueConBlock = ^(NSString * markKey){
```

```
    // Determine connected devices 123456 corresponds to ABCDEF respectively
```

```
    if ([markKey isEqualToString:@"1"]){
```

```
        NSLog(@"Bluetooth device A is connected");
```

```
    }
```

```
    .....
```

```
    else if ([markKey isEqualToString:@"6"]){
```



```

        NSLog(@"Bluetooth device F is connected");
    }
};

// Bluetooth disconnect callback
[Blue4Manager sharedInstance].blueDisBlock = ^(NSString * markKey){
    // Determine disConnected devices
    if ([markKey isEqualToString:@"1"]) {
        NSLog(@"Bluetooth device A is disConnected");
    }

    .....
    else if ([markKey isEqualToString:@"6"]){
        NSLog(@"Bluetooth device F is disConnected");
    }
}

};

//蓝牙数据回调
[Blue4Manager sharedInstance].hzbBlueDataBlock_A = ^(HZLBlueData *blueData, BlueType conBT,
BOOL isFalseCon) {
    if (conBT == BlueType_Pro) {
        if (blueData.bleDataType == BLEMIND) {
            weakSelf.ALabl.text = [NSString stringWithFormat:@"sigal:%d att:%d med:%d
ele:%d ap:%d del:%d theta:%d lowAlp:%d highAlp:%d lowBe:%d highBe:%d lowGa:%d highGa:%d version:%d
grid=%d",blueData.signal,blueData.attention,blueData.meditation,blueData.batteryCapacity,blueData.
ap,blueData.delta,blueData.theta,blueData.lowAlpha,blueData.highAlpha,blueData.lowBeta,blueData.
highBeta,blueData.lowGamma,blueData.highGamma,blueData.hardwareVersion,blueData.grind];
            // when the signal value is 0, the bluetooth device is worn
            // note: if the bluetooth device is connected but not worn, Greater than 0 and less
than or equal to 200
            if(blueData.signal == 0){
                weakSelf.ASignalIV.image = [UIImage
imageNamed:@"signal_zhengChang"];
            }else{
                weakSelf.ASignalIV.image = [UIImage imageNamed:@"signal3.png"];
            }
        }
        else if (blueData.bleDataType == BLEGRAVITY) {
            weakSelf.ACirlLabl.text = [NSString stringWithFormat:@"x:%d y:%d
z:%d",blueData.xvlaue,blueData.yvlaue,blueData.zvlaue];
        }
        else if (blueData.bleDataType == BLERaw) {
            weakSelf.ARawLabl.text = [NSString stringWithFormat:@"raw:%d
eye:%d",blueData.raw,blueData.blinkeye];
        }
    }
};

```

```

    }
}
else if (conBT == BlueType_Jii){
    if (blueData.bleDataType == BLEMIND) {
        weakSelf.ALabl.text = [NSString stringWithFormat:@"sigal:%d att:%d med:%d
ele:%d
ap:%d",blueData.signal,blueData.attention,blueData.meditation,blueData.batteryCapacity,blueData.ap
];

        if(blueData.signal == 0){
            weakSelf.ASignalIV.image = [UIImage
imageNamed:@"signal_zhengChang"];
        }else{
            weakSelf.ASignalIV.image = [UIImage imageNamed:@"signal3.png"];
        }
    }
}
else if (conBT == BlueType_Lite) {
    if (blueData.bleDataType == BLEMIND) {
        weakSelf.ALabl.text = [NSString stringWithFormat:@"sigal:%d att:%d
med:%d del:%d theta:%d lowAlp:%d highAlp:%d lowBe:%d highBe:%d lowGa:%d
highGa:%d",blueData.signal,blueData.attention,blueData.meditation,blueData.delta,blueData.theta,blu
eData.lowAlpha,blueData.highAlpha,blueData.lowBeta,blueData.highBeta,blueData.lowGamma,blueD
ata.highGamma];

        //信号值为 0 即佩戴了蓝牙设备
        //注：如果连接了蓝牙设备而未佩戴，信号值为大于 0 且小于或等于 200
        if(blueData.signal == 0){
            weakSelf.ASignalIV.image = [UIImage
imageNamed:@"signal_zhengChang"];
        }else{
            weakSelf.ASignalIV.image = [UIImage imageNamed:@"signal3.png"];
        }
    }
    else if (blueData.bleDataType == BLERaw) {
        weakSelf.ARawLabl.text = [NSString stringWithFormat:@"raw:%d
eye:%d",blueData.raw,blueData.blinkey];
    }
}

if (isFalseCon) {
    NSLog(@"A device has a false connection");
}

};
.....

[Blue4Manager sharedInstance].hzlblueDataBlock_G = ^(HZLBlueData *blueData, BlueType conBT,
BOOL isFalseCon) {

```

```
.....  
};
```

```
[[Blue4Manager sharedInstance] connectBlue4];
```

```
// Active bluetooth disconnect
```

```
[[Blue4Manager sharedInstance] disconnectBlue4];
```

# iOSBrainLinkProSDKSpecialEdition V1.3

---

## API Reference

### HZLBlueData Reference

#### Overview

The HZLBlueData class is a data model

#### Enum

```
typedef enum : NSUInteger {
    BlueType_NO = 0,
    BlueType_Lite,
    /*The current connection is the BrainLink_Lite data format device with BLEMIND and BLERaw type data*/
    BlueType_Pro,
    /* The current connection is the BrainLink_Pro data format device with BLEMIND, BLEGRAVITY, and BLERaw type data*/
    BlueType_Jii,
    /* The current connection is Jii*/
}BlueType;

typedef NS_ENUM(NSUInteger,BLEDATATYPE){
    BLEMIND    =    0,           // basic brain wave data
    BLEGRAVITY,           // gravity data
    BLERaw,           // blink data
};
```

#### Basic Brainwave Data :

- signal,
- attention,
- meditation,
- delta,
- theta,
- lowAlpha,
- highAlpha,
- lowBeta,
- highBeta,
- lowGamma,
- highGamma,



- ap,
- batteryCapacity,
- hardwareVersion,
- grind

#### Gravity Sensor Data:

- xvlaue,
- yvlaue,
- zvlaue

#### Raw& Blink Data:

- raw,
- blinkeye

#### Note :

When Jii is connected, only signal, attention, meditation, batteryCapacity, ap data type is available.

When BrainLink\_Lite is connected, only signal , attention , meditation , delta , theta , lowAlpha , highAlpha , lowBeta , highBeta , lowGamma , highGamma , raw , blinkeye data type is available.

#### Instructions of some Instance Property

- **signal:** It represents the signal value of the MacroTelligence 's Hardware. When the signal is 0, it means that the MacroTelligence 's Hardware has been put on, and when the signal is greater than 0 and less than or equal to 200, it means that the MacroTelligence 's Hardware is connected to the iPhone.
- **batteryCapacity:** In percentage terms. minimum value is 0, maximum value is 100
- **ap:** Appreciation value
- **hardwareVersion:** Hardware version. The first version value is 255 , when you update the MacroTelligence 's Hardware , the version value will be smaller.
- **xvlaue :** gravity value in The x axis (Pitching Angle)
- **yvlaue :** gravity value in The y axis (Yaw Angle)
- **zvlaue :** gravity value in The z axis (Roll Angle)

## Blue4ManagerReference

### Overview

The Blue4Manager class handles interaction between a MacroTelligence's Hardware and an iOS device.

### Instance Property

#### Successful callback of bluetooth connection

```
@property (nonatomic,copy)Blue4Connect blueConBlock;
```

#### Bluetooth disconnect callback

```
@property (nonatomic,copy) BlueConnectdismiss blueDisBlock;
```

Note: devices are connected in the order of A B C D E F. Six data callbacks (hzluedatablock\_A ...) are used to ensure the independence of data. The data between various devices can be accepted at the same time without mutual influence. Bluetooth 4.0 devices can connect up to six, but it is difficult to connect successfully.

If you want to use a single connection, the input parameter of ableddevicesum is 1. Only call hzluedatablock\_A.

#### Data callback for every device.

```
@property(nonatomic,copy)Blue4DataBlock hzlblueDataBlock_A;  
@property(nonatomic,copy)Blue4DataBlock hzlblueDataBlock_B;  
@property(nonatomic,copy)Blue4DataBlock hzlblueDataBlock_C;  
@property(nonatomic,copy)Blue4DataBlock hzlblueDataBlock_D;  
@property(nonatomic,copy)Blue4DataBlock hzlblueDataBlock_E;  
@property(nonatomic,copy)Blue4DataBlock hzlblueDataBlock_F;
```

#### Connection status of every device

```
@property (nonatomic,assign)BOOL connected_A;  
@property (nonatomic,assign)BOOL connected_B;  
@property (nonatomic,assign)BOOL connected_C;  
@property (nonatomic,assign)BOOL connected_D;  
@property (nonatomic,assign)BOOL connected_E;  
@property (nonatomic,assign)BOOL connected_F;
```

### Method

#### Print log does not print by default

```
+ (void)logEnable:(BOOL)enable;
```

#### Initialization (singleton)

```
+ (instancetype)shareInstance;
```

### Parameter Configuration :

#### Parameter interpretation :

blueNames: **Device name (Bluetooth 4.0 device) or MFI (Bluetooth 3.0 device) that can be connected**

```
NSArray *blueNames = @[@"BrainLink",@"BrainLink_Pro",@"jii@jii-***"];
```

1.jii@jii-\*\*\* : indicates the device name with the prefix of jii - that can be connected , 'jii@' means jii device. The content After @ is the device name. '\*\*\*' indicates the prefix is the same.

ableDeviceSum: Number of Bluetooth devices that can be connected

```
-(void)configureBlueNames:(NSArray *)blueNames ableDeviceSum:(int)deviceSum;
```

#### Connect bluetooth device

```
-(void)connectBlue4;
```

#### Disconnect bluetooth device

```
-(void)disConnectBlue4;
```

**Manual test of false connection (definition of false connection: when signal is equal to 0 and the continuous 10 values of attention and mediation remain unchanged, it is considered to be a false connection. SDK will disconnect the Bluetooth connection of the current device and automatically connect again)**

```
-(void)testAFalseCon:(BOOL)isTest; //Manual test A device false connection
```

```
-(void)setTestToZero;//Cancel all manual test false connections
```